



-  **Video output port**
-  **Micro SD card slot**
Supporting max 32G
-  **Mini USB port**
For charging



Operation board
Providing for snap shot and video record function.

TRACHWAY[®]

VIDEO INTUBATION SERIES

One Handle but TWO+ Choices



Wheel button
Use one finger to adjust the angle of monitor, 60° upward and 45° downward.



Dis assemble latch
In order that the monitor can be easily slide-out from handle, lift up the latch



Main Power



Quick connector
Fits with Stylet series and Laryngoscope series



Built-in wireless transmitter.
Transmitter and receiver will activate simultaneously when the device been separated while power in turned on status.



Line with conventional



Water-proof




Immerse Disinfection



Foggy free

Video Laryngoscope

- Same as conventional Miller and Mac. Blade
- Same skill but non direct view to the vocal cord
- No need special curve bending on ETT
- Upgrade the C&L at least one grade to reduce incidence of difficult airway.
- Foggy Free
- No consumption

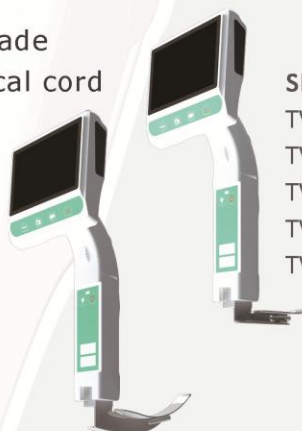

SIZES-

- TVL-4220 Miller #0
- TVL-4221 Miller #1
- TVL-4232 Macintosh #2
- TVL-4233 Macintosh #3
- TVL-4234 Macintosh #4

Video Light Stylet



Red LED on Tip



Malleable



Water-proof



Immerse Disinfection




Foggy free



- Malleable stylet can fit difficult airway patients.
- Limited mouth opening and neck movement
- Atraumatic tip without concave for patient safety.
- Cervical spine immobilization intubation
- No sore throat
- Foggy free
- Water proof for high level disinfection
- Good for General ETT,
Preform ETT,
2 Lumen ETT intubation

SIZES-

TVI-4102	5.0	345mm	Regular
TVI-4105	5.0	420mm	Double Lumen
TVI-4106	4.0	300mm	Pediatric
TVI-4108	4.0	420mm	Double Lumen

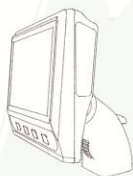

TVI-4050
3" TFT Monitor
W 32 x H 110 x D 39 (mm)
110 g



TVI-4051
Rechargeble Handle
W 90 x H 110 x D 88 (mm)
160g